

Datasheet for 600-143-096

Fluorescein Antibody DyLight™ 649 Conjugated

Overview

Description:	Anti-Fluorescein (GOAT) Antibody DyLight™ 649 Conjugated - 600-143-096
Item No.:	600-143-096
Size:	100 µg
Applications:	Dot Blot, WB
Reactivity:	Fluorescein
Host Species:	Goat

Product Details

Background:	Fluorescence technology is widely used to detect proteins. However, many common visible fluorophores often result in considerable background fluorescence in the visible range. Visible fluorophores are rarely used for membrane-based protein detection because of this high background. The emission spectra for this DyLight™ conjugate match the principle output wavelengths of most common fluorescence instrumentation. DyLight™ antibody and reagent conjugates are specifically designed for protein detection methods that cover much of the visible spectrum and into the near-infrared region to visualize proteins in western blotting and other applications. Very low background fluorescence in the IR range provides for a much higher signal-to-noise ratio than visible fluorophores. Detection levels in the picogram range on Western blots rival the sensitivity of chemiluminescence on film. DyLight™ conjugates provide an ultra-sensitive and convenient alternative to standard chemiluminescent protein detection methods, as well as a valuable tool for multicolor imaging.
Synonyms:	goat Anti-FITC antibody DyLight™ 649 conjugation, DyLight™ 649 conjugated goat anti-fluorescein antibody, DyLight™ 649 conjugated goat Anti-FITC antibody, goat anti-fluorescein antibody DyLight™ 649 conjugation
Host Species:	Goat
Conjugate:	DyLight™ 649
Clonality:	Polyclonal
Format:	IgG
F/P Ratio:	3.7

Target Details

Reactivity:	Fluorescein
Immunogen Type:	Other
Immunogen:	Fluorescein conjugated to Goat IgG
Purity/Specificity:	Fluorescein Antibody Dylight™ 649 Conjugated was prepared from monospecific antiserum by immunoaffinity chromatography using Fluorescein IgG coupled to agarose beads followed by solid phase adsorption(s) to remove any unwanted reactivities. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Goat Serum and Fluorescein conjugated Bovine Serum Albumin.

Application Details

Tested Applications:	Dot Blot, WB
Application Note:	Fluorescein Antibody Dylight™ 649 Conjugated has been tested by dot blot and western blot and is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA) and fluorescent western blotting. This product is also suitable for multiplex analysis, including multicolor imaging, utilizing various commercial platforms.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
FLISA:	>1:10,000
IF:	>1:5,000

Formulation

Physical State:	Lyophilized
Concentration:	1.0 mg/mL by UV absorbance at 280 nm
Buffer:	0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2
Preservative:	0.01% (w/v) Sodium Azide
Stabilizer:	10 mg/mL Bovine Serum Albumin (BSA) - Immunoglobulin and Protease free
Reconstitution Volume:	100 µL
Reconstitution Buffer:	Restore with deionized water (or equivalent)







Shipping & Handling

Shipping Condition:	Ambient
Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.
Expiration:	Expiration date is one (1) year from date of receipt.

Images


Diagram


Properties of DyLight™ Conjugates.

Emission	Color	DyLight™ Dye	Ex/Em (nm)	ϵ (M ⁻¹ cm ⁻¹)	Similar Dyes
Blue		405	400/420	30,000	Alexa™ 405, Cascade Blue
Green		488	493/518	70,000	Alexa™ 488, Cy2®, FITC
Yellow		549	550/568	150,000	Alexa™ 546, Alexa 555, Cy3®, TRITC
Red		649	646/674	250,000	Alexa™ 647, Cy5®
Near Infrared		680	682/715	140,000	Alexa™ 680, Cy5.5®, IRDye™ 700
Infrared		800	770/794	270,000	IRDye™ 800

Western Blot

DyLight™ dyes can be used for two-color Western Blot detection with low background and high signal. Anti-tubulin was detected using a DyLight™ 549 conjugate. Anti-TNF α was detected using a DyLight™ 649 conjugate. The image was captured using the Typhoon™ 9410 Imaging System.

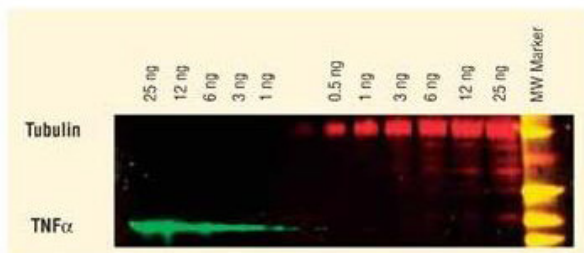

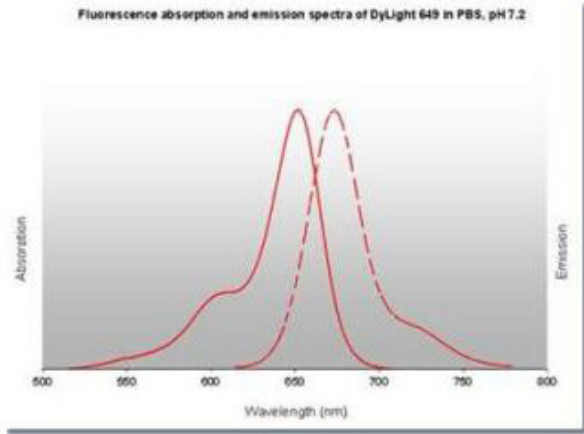**Diagram**

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.